I'm a genuine  
La-Z-Boy® product

---

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---

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Your opinions are very important to us. Complete our online **survey** within 30 days at **www.lzb-delivery.com** and your name will be entered into a drawing for \$2500 worth of La-Z-Boy furniture. No purchase necessary to win.

Void in Quebec and where prohibited by law. Odds of winning depend on number of entries received.

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L A Z B O Y®

## Important Safety Instructions:

Important: Read these instructions thoroughly before using.

### **DANGER**

– To reduce the risk of electric shock, serious injury or death:

- Do not let children play on or near the chair when in use. Serious injury or death could occur if head or body parts are under the chair arms when the chair is moved from the full lift to the seated position.
- Disconnect the electrical plug from outlet before cleaning, moving, servicing or when not in use for an extended period of time.
- Do not use this product if it is wet, in water, around water or any liquid. Do not handle the plug or remote control with wet hands.
- To prevent suffocation keep plastic bags away from children.
- Plug directly into a wall outlet. Do not use an extension cord, adapter or surge protector.

### **WARNING**

– To reduce the risk of electric shock or serious injury:

- Sit down in the unit slowly. To avoid tip over, do not throw your weight against the back of the unit.
- Keep hands, fingers and body parts away from the mechanism. To reduce the risk of head entrapment and other injuries, keep children and pets away from openings and moving parts in and around furniture.
- Only the occupant should operate the unit. Never attempt to open the legrest by pulling on it.
- Do not stand, sit, or place extra weight on the arms, back or open legrest.
- To safely exit the unit, bring the back upright and close the legrest.
- Do not use this product if it is damaged. Do not use attachments or substitute parts.
- Check the condition of the power transformers, electrical cords, plugs, and connectors monthly. Keep the power transformers and electrical cords away from heat sources. Never operate the unit with a damaged power transformer, electrical cord, plug, or connector. Contact your La-Z-Boy® dealer for service.
- Do not use unit if power control is damaged or wet. Do not use unit if the hand wand has been dropped in water or any liquid. Contact your La-Z-Boy® dealer for service.
- Some medications and physical limitations may impair your ability to operate this product in a safe manner. Consult your physician if you are taking prescribed or over the counter medication, or if you have physical limitations.

- Users of medical devices should contact their health care provider before using this product. Do not use this product if you have poor circulation, are immobile, or cannot sense heat, except as directed by your health care provider.
- Do not use this product where aerosol (spray) products are being used or where oxygen is being administered.
- In the event of a power failure, the battery backup provided may not lift the user to a fully raised position. To ensure lift functionality during a power outage, a battery backup system must be purchased separately. Contact your La-Z-Boy® dealer for more information. Replace the power transformer batteries every six months and after each power failure use.
- Before moving the product, disconnect the electrical plug from the wall outlet and secure the mechanism to prevent unexpected opening or movement. Never move the product by the mechanism, legrest, handle, reclining back or swivel base.
- Do not attempt to disassemble or service the pump or motor; no serviceable parts are inside. Contact your La-Z-Boy® dealer for service.
- The 1HL/1ML lift chairs are intended for a single occupant, not exceeding 500 pounds. Do not operate the chair with more than one occupant in or on any part of the chair.
- The 1LF lift chairs are intended for a single occupant, not exceeding 300 pounds. Do not operate the chair with more than one occupant in or on any part of the chair.
- Do not allow electrical cord to be pinched between moving parts or between the frame and the floor. Do not place cord under carpet or in location where it could become worn.
- Do not use pillow, blanket, or other cover between occupant and chair surface when heat is engaged.
- Burns may result from improper use and care of internal heating pads. To avoid risk of burns, do not use pins, do not allow heating pad area to become wet, and do not remove upholstery cover.

### **NOTICE**

– To reduce the risk of product damage:

- To prevent damage to the power transformer, do not plug in until instructed.

Save These Instructions

# Assembly Instructions:

## **WARNING**

– To reduce the risk of injury:

- Do not operate the unit until the back is installed and locked.

### Assembly Instructions

1. Carefully cut the cable tie securing the lift base (FIGURE 1).
2. Align the brackets on the back with the brackets on the body (FIGURE 2).
  - The back brackets may be located on the outside or inside of the back assembly.
3. Lower the back onto the body, making sure the back and body brackets are completely engaged.
4. Push the locking levers down to lock the brackets (FIGURE 3).
  - If you are unable to install the back properly, contact your La-Z-Boy dealer for assistance.

FIGURE 1

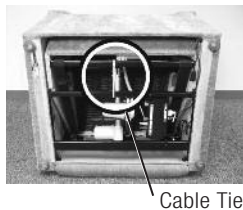
FIGURE 2

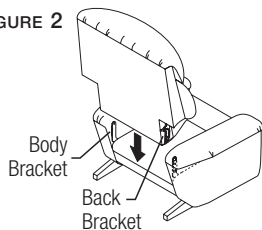
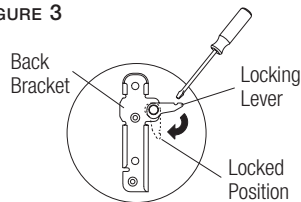
FIGURE 3

5. Connect the back massage connectors for 1ML styles, and the back heating pad connectors for the 1HL and 1ML styles (FIGURE 4).
6. Connect the Power Transformer connectors (FIGURE 5).
  - Remove the Power Transformer from the box. Connect the Power Transformer to the connector behind the chair. Do not plug in until instructed.
7. Attach the black non-woven flap to form a pocket between the seat and the back.
  - Push the black non-woven flap between the seat and the back, toward the back of the chair. Roll chair forward to rest on front arms.
  - Attach the black non-woven flap to the back seat spring rail by pressing the fastening tape together. Roll chair back on to base.
8. Plug the power transformer directly into an outlet.
  - Raise the lift chair. Attach the fastening tape on the back upholstery cover to the fastening tape on the inner body back panel.
9. Proceed to the online Operating Instructions.

FIGURE 4

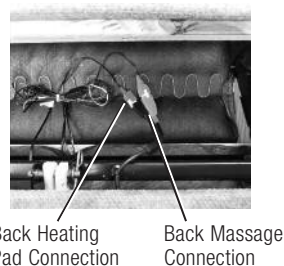
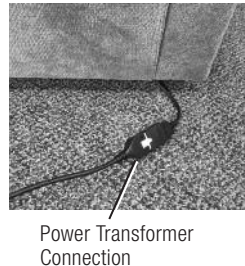
FIGURE 5

Proceed to the Operating Instructions at [la-z-boy.com/Support/Product-Manuals](http://la-z-boy.com/Support/Product-Manuals) for complete assembly instructions, operating instructions, and safety instructions.

For fabric and leather care and cleaning instructions, visit [la-z-boy.com/Support/Fabric-Care](http://la-z-boy.com/Support/Fabric-Care)

If you do not have internet access, contact your La-Z-Boy dealer. If you are unable to reach your La-Z-Boy dealer, call La-Z-Boy Comfort Care at 855-80-COMFORT (or 855-802-6636).

## Warranty Information:

La-Z-Boy provides the following warranties, subject to the stated limitations.

### LIMITED WARRANTY – FABRICS

Fabrics are warranted for one (1) year against seam failure, fabric separation and nap loss. Pilling, dye lot match, and fading are not covered. Labor: One (1) year.

For directions on care and cleaning of fabric, visit: [la-z-boy.com/Support/Fabric-Care](http://la-z-boy.com/Support/Fabric-Care) or call 855-80-COMFORT (or 855-802-6636). Improper cleaning/treatment or application of treatment voids your warranty.

### LIMITED WARRANTY – SEAT CUSHIONS

Standard foam cushions have a Limited Warranty for one (1) year.

Labor for attached seat cushions is one (1) year.

### LIMITED WARRANTY – PARTS

Certain parts are warranted against manufacturing and material defects, as shown below.

Labor: One (1) year unless noted. What is covered:

- Frame: Lifetime of unit
- Base: Three (3) years
- Back and Seat Springs: Three (3) years
- Recliner/Lift mechanism: Three (3) years

### LIMITED WARRANTY – ELECTRICAL SYSTEMS

Electric motor and electrical systems have a three (3) year Limited Warranty.

Labor: One (1) year.

### LIMITATIONS APPLICABLE TO ALL WARRANTIES

ALL WARRANTIES ARE TO ORIGINAL PURCHASER FROM AUTHORIZED DEALERS ONLY. WARRANTIES ARE NOT TRANSFERABLE.

ALL WARRANTIES ARE LIMITED TO MANUFACTURING AND MATERIAL DEFECTS THAT APPEAR UNDER NORMAL INDOOR RESIDENTIAL USE. ADDITIONAL LIMITATIONS DESCRIBED HEREIN APPLY. WARRANTIES ARE FOR REPAIR, REPLACEMENT, OR SUBSTITUTION ONLY, IN LA-Z-BOY'S SOLE DISCRETION. REFUNDS ARE NOT AVAILABLE.

LA-Z-BOY WILL PAY REASONABLE AND CUSTOMARY RATES FOR LABOR DURING THE PERIOD OF THE LABOR WARRANTY. YOU WILL BE RESPONSIBLE FOR LABOR COSTS AFTER THE STATED LABOR TIME PERIOD.

YOU ARE RESPONSIBLE FOR IN-HOME INSPECTION FEES, TRANSPORTATION OF PRODUCT TO AND FROM YOUR DEALER AND TO AND FROM THE LA-Z-BOY FACTORY. (MOST REPAIRS DO NOT REQUIRE RETURN TO THE FACTORY.)

WARRANTIES ARE VOIDED BY EVIDENCE OF EXCESSIVE SOILING, IMPROPER CLEANING OR TREATMENT, ABUSE, OR ABNORMAL USE. WARRANTIES DO NOT APPLY TO PRODUCTS USED FOR RENTAL, BUSINESS, COMMERCIAL, INSTITUTIONAL, OR OTHER NON-RESIDENTIAL USES.

**THESE WARRANTIES ARE YOUR EXCLUSIVE REMEDY AND NO OTHER WARRANTIES, EXPRESS OR IMPLIED, INCLUDING, BUT NOT LIMITED TO, THE IMPLIED WARRANTIES OF MERCHANTABILITY AND FITNESS FOR A PARTICULAR PURPOSE, APPLY. IN NO EVENT WILL LA-Z-BOY BE LIABLE FOR INCIDENTAL OR CONSEQUENTIAL DAMAGES, EVEN IF IT HAD REASON TO KNOW IN ADVANCE THAT SUCH DAMAGES WERE POSSIBLE. THESE WARRANTIES GIVE YOU SPECIFIC LEGAL RIGHTS, AND YOU MAY HAVE OTHER RIGHTS THAT VARY FROM JURISDICTION TO JURISDICTION. SOME JURISDICTIONS DO NOT ALLOW THE EXCLUSION OF IMPLIED WARRANTIES OR THE LIMITATION OF DAMAGES, SO THE ABOVE EXCLUSIONS AND LIMITATIONS MAY NOT APPLY TO YOU. IN THE EVENT THE WARRANTY EXCLUSIONS DO NOT APPLY TO YOU, THEN, WHERE LEGALLY PERMISSIBLE, SUCH WARRANTIES ARE LIMITED IN DURATION TO THE WARRANTY PERIOD AND NO WARRANTIES APPLY AFTER THAT PERIOD.**

### CLAIMS

To make a claim under this warranty, please contact the authorized dealer from whom you made your purchase. (If your original dealer is not available, please contact another authorized dealer with your proof of purchase.) La-Z-Boy products may not be altered by any person, dealer, or company without the express written authorization of La-Z-Boy.

If you feel, after presenting your warranty claim to your dealer, that it has not been resolved, you may write to:

La-Z-Boy Incorporated  
1284 N. Telegraph Road  
Monroe, MI 48162  
Attn: Comfort Care

Or visit, [www.la-z-boy.com/ContactUs](http://www.la-z-boy.com/ContactUs)

Please include the dealer's name and address, furniture delivery date, ALI # (Acknowledgment-Line-Item #), style, cover number, and a brief explanation of the problem.

**PRODUCT REGISTRATION**

To register your product visit:  
**[la-z-boy.com/Support/Product-Registration](http://la-z-boy.com/Support/Product-Registration)**

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©2014 La-Z-Boy Incorporated